

BLIMS

Bounty Laboratory Information Management System

User Manual

Introduction

BLIMS (Bounty Laboratory Information Management System) is a web-based system designed to help laboratory personnel record, manage, and monitor laboratory test results efficiently. This manual is intended for non-technical users responsible for recording laboratory sample tests.

System Objective.

To modernize request handling, testing workflows, and data security by replacing manual, paper-based processes with a streamlined, digital platform.

User Roles & Access

Laboratory users can log in, view assigned samples, record and update test results, and review test history. Access depends on permissions assigned by the system administrator.

System Requirements

- Desktop or laptop computer
- Internet connection (with VPN)
- Google Chrome (recommended)

Logging In to BLIMS

1. Open your web browser.
2. Enter the BLIMS system URL.
3. Enter your username and password.
4. Click Login.

Registration

Transactional
Registration

Add Registration

All Fields
Search here

Job Order No: PFMQA-JN-0126-0001
Add Detail
Samples 2

No.	Laboratory Code	Sample Name	Production ...	Shipment Supplier	Plate / ...	Batch / ...	Type o...	Labora...	Test P...	Test Status	Reason	COA R...	Remarks	Action
1	29-0126-0001	ADISODIUM	2026-01-02	A & S INTERNATIONAL	1221	2001	FEEDS AND RAW MATERIALS	MOISTURE CONTENT	MOISTURE, %	ONGOING	No Reason Indicated	No	No Remarks	
2	29-0126-0002	ADISODIUM	2026-01-02	A & S INTERNATIONAL	ABB 5515	2001	FEEDS AND RAW MATERIALS	CRUDE PROTEIN	CRUDE PROTEIN, %	ONGOING	No Reason Indicated	Yes	No Remarks	

Click “Add Registration” to create a new sample registration then you will be redirect to creating registration module.

Registration

Transactional
Registration

Registration

Internal Feedmill *

Select

Commercial Feedmill

Select Feed

Address

Enter address

External Client

Contact Number *

0917 123 4567

Nutritionist *

Select

Laboratory Location *

Select Location

Delivery Type

Courier

Internal Handling

Way Bill Number

Enter Way Bill Number

Submitted By

JAMES DIMACULANGAN

Date Submitted

2026-01-02

Remarks

Enter any remarks

Attach File

Choose File

No file chosen

Sample Details *

Sample Name	Type of Sample	Test Code	Production/Delivery D...	Shipment/Supplier	Plate/Van Number	Batch/Lot Number	Type of Lead Time	COA	Action
Sample	Type of Sample	Test Code	dd/mm/yyyy	Supplier	Plate/Van	Batch/Lot Number	Select Lead Time		+ -

Save

Cancel

Here just fill up the all the dropdown details required for all fields.

Note. To enabled the Sample Details, you need to select the feedmills, nutritionist and Laboratory Location.

Sample Details *


Sample Name	Type of Sample	Test Code	Production/Delivery D...	Shipment/Supplier	Plate/Van Number	Batch/Lot Number	Type of Lead Time	COA	Action
Sample	Type of Sample	Test Code	dd/mm/yyyy	Supplier	Plate/Van	Batch/Lot Number	Select Lead Time		+ -

Save

Cancel

Each field in the sample details is a dropdown fields, just select all the fields and also, we

have an action button where we can add other samples for this Registration. Click “Save” to submit. Once save Success Message will appear.

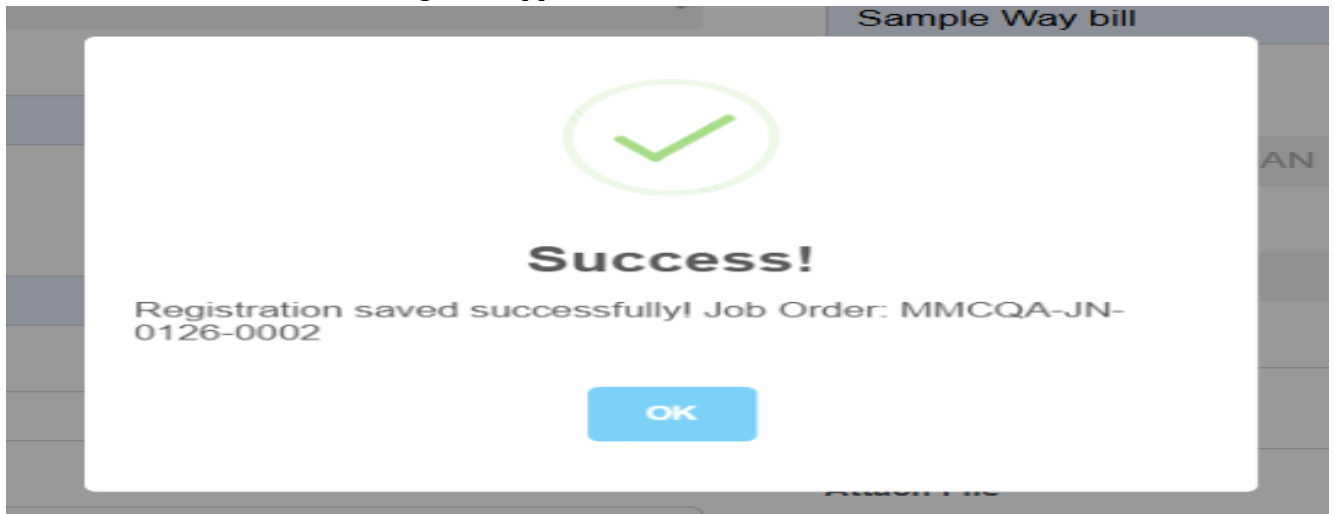

Then go back to View of Created Registration. Here we can edit a registered samples(*only if you have an access to do so*).

Job Order No: MMCQA-JN-0126-0002

+ Add Detail

Samples: 1

No.	Laboratory Code	Sample Name	Production ...	Shipment Supplier	Plate / ...	Batch / ...	Type o...	Labora...	Test P...	Test Status	Reason	COA R...	Remarks	Action
1	12-0126-0003	ARGENTINE CORN	2026-01-02	A & S INTERNATIONAL	AJA 9142	B#C2523 -07	FEEDS AND RAW MATERIALS	CRUDE FAT	CRUDE FAT, %	ONGOING	No Reason Indicated	Yes	No Remarks	

Click “Add Detail” if you want to add a new samples to existing Job Order then modal will appear fill up all the details then Click “Save” if done.

Add Sample Detail

Sample Name

Select

Type of Sample

Select

Test Code

Select

Production/Delivery Date

dd/mm/yyyy

Shipment/Supplier

Select

Plate / Van Number

Select

Batch / Lot Number

Select

Type of Lead Time

COA Required:

☐

Save

Cancel

If you want to append a new sample with sample lab code under the same Job Order, Click the Action button “replicate icon” then modal will appear and it will displayed the same details you want to replicate.

Job Order No: MMCGA-JN-0126-0002

+ Add Detail

Samples: 1

No.	Laboratory Code	Sample Name	Production ...	Shipment Supplier	Plate / ...	Batch / ...	Type o...	Labora...	Test P...	Test Status	Reason	COA R...	Remarks	Action
1	12-0126-0003	ARGENTINE CORN	2026-01-02	A & S INTERNATIONAL	AJA 9142	B#C2523 -07	FEEDS AND RAW MATERI ALS	CRUDE FAT	CRUDE FAT, %	ONGOING	No Reason Indicated	Yes	No Remarks	

Replicate Sample Details: 12-0126-0003

Sample Name

ARGENTINE CORN

Type of Sample

Select

Test Code

Select

Production/Delivery Date

02/01/2026

Shipment/Supplier

A & S INTERNATIONAL

Plate / Van Number

AJA 9142

Batch / Lot Number

B#C2523-07

Type of Lead Time

COA Required:

☒

Save

Cancel

Then fill up all the required fields then Click “Save”, then modal for confirmation will appear, click “Yes, Proceed”.

Replicate Sample Details: 12-0126-0003

Sample Name

ARGENTINE CORN

Type of Sample

LIMESTONE

Test Code

NIR-CF

Production/Delivery Date

02/01/2026

Shipment/Supplier

A & S INTERNA

Plate / Van Number

AJA 9142

Batch / Lot Number

B#C2523-07

Type of Lead Time

Regular

COA Required:

☒

Save

Cancel

Replicating this Sample


Are you sure you want to Replicate this Samples?

✓ Yes, Proceed

✕ Cancel

Verification Module

Allows users to view and update the verification status of registered samples, categorized as Ongoing, Cancelled, or Failed. Users can assign verification results and provide reasons or remarks when required.

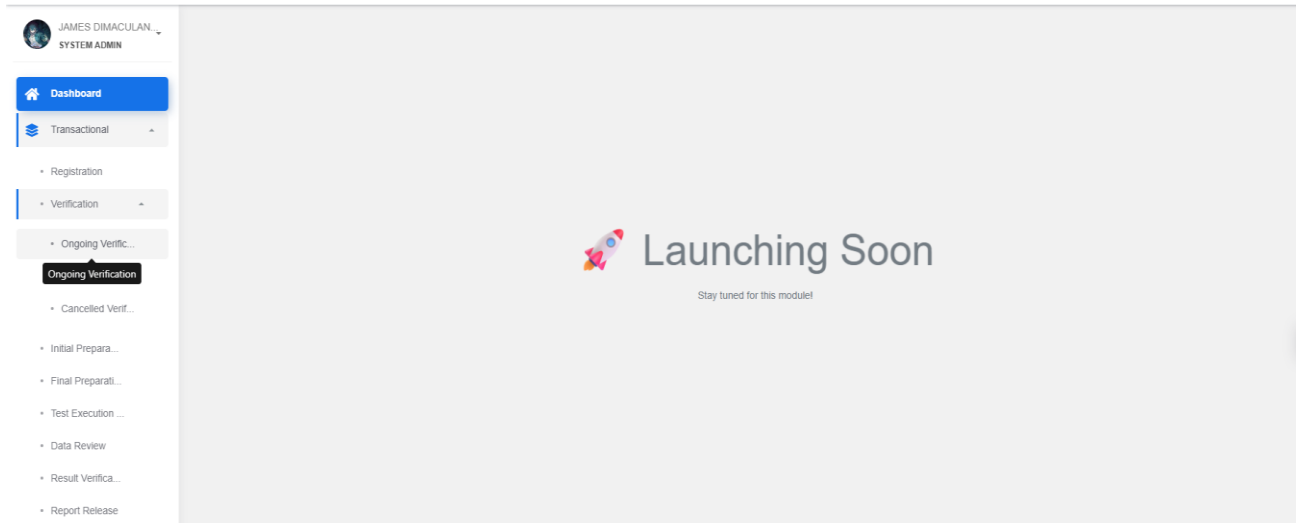

Click Ongoing Verification Module.

Ongoing Verification

Job Order No: MMCQA-JN-0126-0002

Samples: 2

No.	Laboratory Code	Sample Name	Production ...	Shipment Supplier	Plate / V...	Batch / L...	Type of ...	Laborato...	Test Par...	Test Status	Reason	COA Req...	Remarks
1	12-0126-0003	ARGENTINE CORN	2026-01-02	A & S INTERNATIONAL	AJA 9142	B#C2523-07	LIMESTON E	NIR PROXIMAT E	CRUDE FIBER, %	Select	Select	Yes	
2	12-0126-0003	ARGENTINE CORN	2026-01-02	A & S INTERNATIONAL	AJA 9142	B#C2523-07	FEEDS AND RAW MATERIAL S	CRUDE FAT	CRUDE FAT, %	Select	Select	Yes	

Select test status and reasons for this. *If On hold and failed the selected test status, reasons and remarks field are required. It will also notify the User who created the Registration via email.* Then Click "Save" then Confirmation modal will appear.

Confirm Verification

Status Summary

Passed
1

On Hold
0

Failed
0

Are you sure you want to proceed with verification?

Yes, Proceed

Cancel

Click “Yes, proceed”.

Failed Verification

All Fields
Search here

Job Order No: PFMQA-JN-1125-0054

Samples 1

No.	Laboratory Code	Sample Name	Production ...	Shipment Supplier	Plate / V...	Batch / L...	Type of ...	Laborato...	Test Par...	Test Status	Reason	COA Req...	Remarks
1	16-1125-0068	ARGENTINE SOYA	2025-11-26	AAFCO	ABB 5515	B#C2523-08	MOLASSE S	NIR PROXIMATE	ASH, %	FAILED	DAMAGED OR IMPROPER PACKAGING	No	failed

Note: If samples are not edit or change test status after 3 business days it will automatically Cancelled by our system.

Cancelled Verification

All Fields
Search here

Job Order No: MMCQA-JN-1125-0099

Samples 3

No.	Laboratory Code	Sample Name	Production ...	Shipment Supplier	Plate / V...	Batch / L...	Type of ...	Laborato...	Test Par...	Test Status	Reason	COA Req...	Remarks
1	16-1125-0085	ADISODIUM	2025-11-27	A & S INTERNATIONAL	1221	2001	LIMESTONE	NIR PROXIMATE	CRUDE FIBER, %	CANCELLED	DAMAGED OR IMPROPER PACKAGING	No	Auto-cancelled by system after 3 business days
2	16-1125-0089	ADISODIUM	2025-11-27	A & S INTERNATIONAL	1221	B#L2323-02	LIMESTONE	NIR PROXIMATE	CRUDE FIBER, %	CANCELLED	INCORRECT SAMPLE WEIGHT	No	Auto-cancelled by system after 3 business days
3	16-1225-0002	ADISODIUM	2025-11-27	ACX DIAMOND	1221	B#C2523-07	OIL	CA(EDTA)	CALCIUM, %	CANCELLED	INSUFFICIENT RESOURCES (E.G., CHEMICAL REAGENTS, CONSUMABLES, ETC.)	No	Auto-cancelled by system after 3 business days

Initial Preparation Module

The Initial Preparation module allows users to manage and edit the pre-analytical process before moving samples to Final Preparation.

Initial Preparation Transactional > Initial Preparation

All Fields Search here

Job Order No: MMCQA-JN-0126-0002 Samples: 1

No.	Laboratory Code	Sample Name	Production ...	Shipment Supplier	Plate / Van Number	Batch / Lot Number	Type of Sample	Laboratory Tests	Test Parameter	Pre Analytical Process
1	12-0126-0003	ARGENTINE CORN	2026-01-02	A & S INTERNATIONAL	AJA 9142	B#C2523-07	LIMESTONE	NIR PROXIMATE	CRUDE FIBER, %	Select

Select Pre-analytical Process for this Sample then Click “Save” then confirmation modal will appear. Click “Yes, proceed”

Confirm Initial Preparation

Pre Analytical Process

Homogenizing

1

Grinding

0

As Is

0

Are you sure you want to proceed with Initial Preparation?

✓ Yes, Proceed **✗ Cancel**

Final Preparation Module

The Final Preparation module allows users to manage and edit the preparation Verification before moving samples to Test Execution and Data Entry.


Final Preparation Transactional > Initial Preparation > Final Preparation

All Fields Search here

Job Order No: MMCQA-JN-0126-0002 Samples: 1

No.	Laboratory Code	Sample Name	Production ...	Shipment Suppl...	Plate / Van Num...	Batch / Lot Nu...	Type of Sample	Laboratory Tests	Test Parameter	Pre Analytical A...	Preparation Ver...	Remarks
1	12-0126-0003	ARGENTINE CORN	2026-01-02	A & S INTERNATIONAL	AJA 9142	B#C2523-07	LIMESTONE	NIR PROXIMATE	CRUDE FIBER, %	HOMOGENIZING	Select	

Select in the dropdown field. *If failed remarks are required and it will notify the user via email and the samples will proceed to Initial Preparation.* Click “Save” then formation modal will appear.

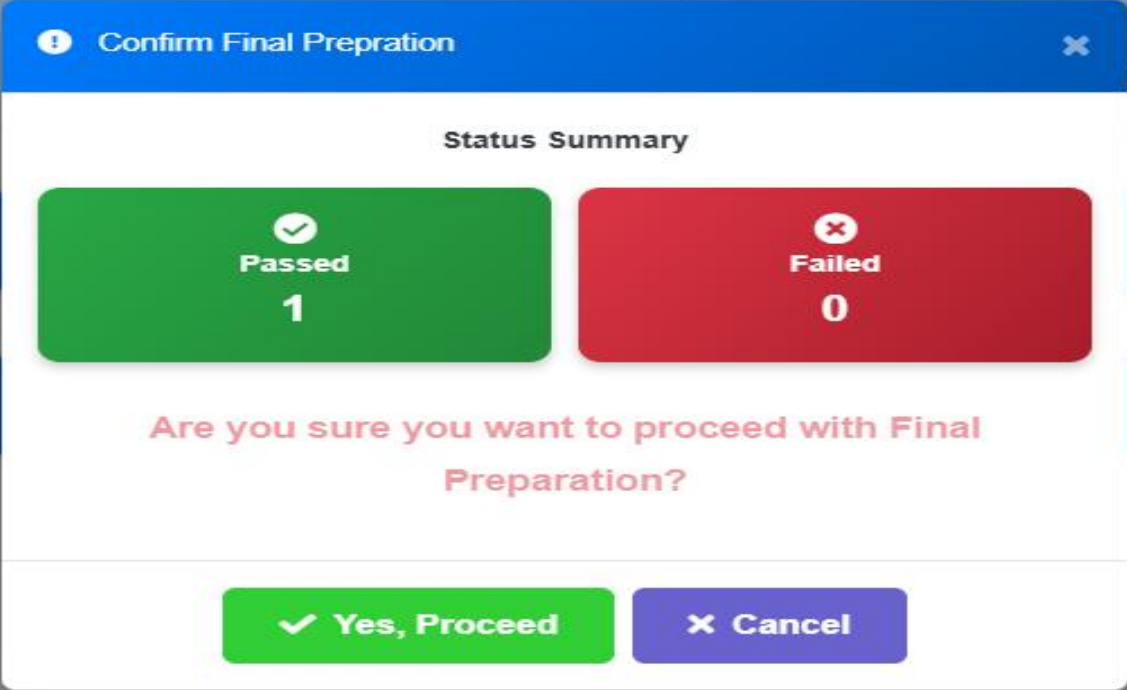


The image shows a modal window titled "Confirm Final Preparation". It features a "Status Summary" section with two boxes: a green box labeled "Passed" with a checkmark icon and the number "1", and a red box labeled "Failed" with an 'X' icon and the number "0". Below this, a question asks, "Are you sure you want to proceed with Final Preparation?". At the bottom, there are two buttons: a green "Yes, Proceed" button and a purple "Cancel" button.

Click “Yes, proceed” to move the samples to Test Execution and Data Entry.

Test Execution and Data Entry Module

The Test Execution and Data Entry module allows users to manage and edit the status also provide the Lab Result and Remarks; it also displayed here the Lead time.

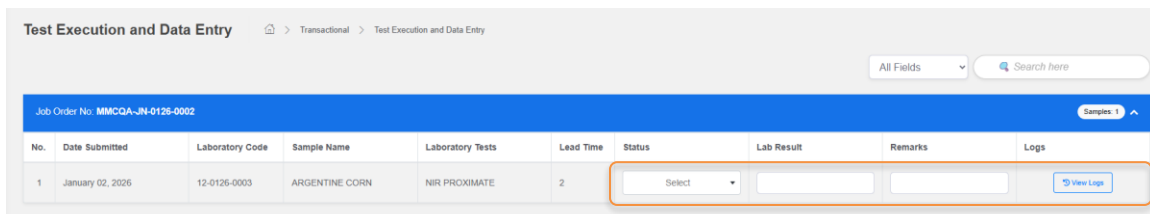

The image shows the "Test Execution and Data Entry" interface. It includes a breadcrumb trail: "Test Execution and Data Entry > Transactional > Test Execution and Data Entry". There is a search bar with "All Fields" and "Search here". Below this is a table with the following data:


No.	Date Submitted	Laboratory Code	Sample Name	Laboratory Tests	Lead Time	Status	Lab Result	Remarks	Logs
1	January 02, 2026	12-0126-0003	ARGENTINE CORN	NIR PROXIMATE	2	Select			View Logs

Select for dropdowns.

If the selected is Completed in the status Lab result is required.

If the Sample Beyond Lead time it will notify Manager, Supervising Chemist and Analyst via email and Remarks is required.

Click “View logs” to view samples history or record.

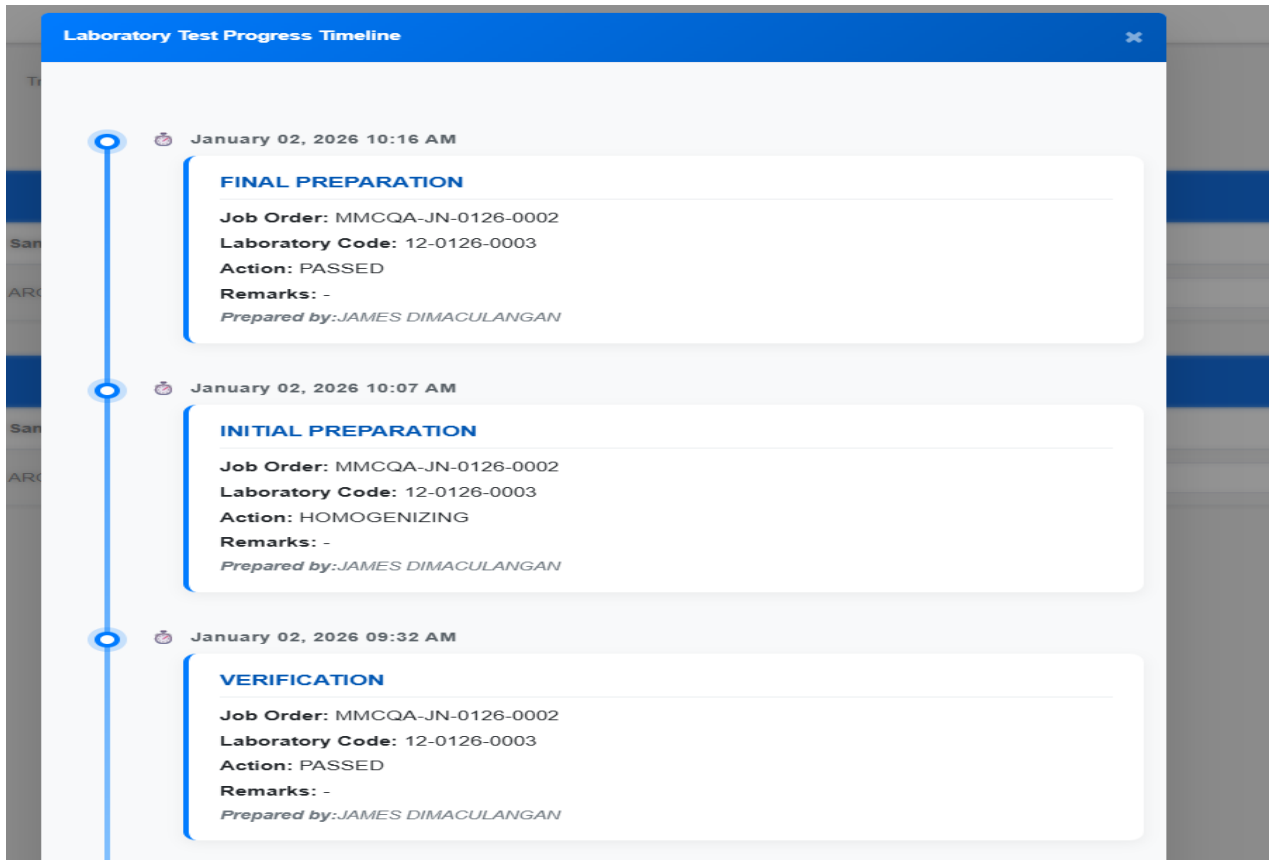

The screenshot displays a 'Laboratory Test Progress Timeline' window. It features a vertical timeline on the left with three circular markers. Each marker corresponds to a test event. The first event, dated January 02, 2026 at 10:16 AM, is titled 'FINAL PREPARATION'. The second event, dated January 02, 2026 at 10:07 AM, is titled 'INITIAL PREPARATION'. The third event, dated January 02, 2026 at 09:32 AM, is titled 'VERIFICATION'. Each event card contains the following details: Job Order: MMCQA-JN-0126-0002, Laboratory Code: 12-0126-0003, Action: PASSED, Remarks: -, and Prepared by: JAMES DIMACULANGAN.

Laboratory Test Progress Timeline

January 02, 2026 10:16 AM

FINAL PREPARATION

Job Order: MMCQA-JN-0126-0002
Laboratory Code: 12-0126-0003
Action: PASSED
Remarks: -
Prepared by: JAMES DIMACULANGAN

January 02, 2026 10:07 AM

INITIAL PREPARATION


Job Order: MMCQA-JN-0126-0002
Laboratory Code: 12-0126-0003
Action: HOMOGENIZING
Remarks: -
Prepared by: JAMES DIMACULANGAN

January 02, 2026 09:32 AM

VERIFICATION

Job Order: MMCQA-JN-0126-0002
Laboratory Code: 12-0126-0003
Action: PASSED
Remarks: -
Prepared by: JAMES DIMACULANGAN

Click “Save”, then confirmation modal will appear click “Yes, Proceed” to move the samples for Data Review Module.

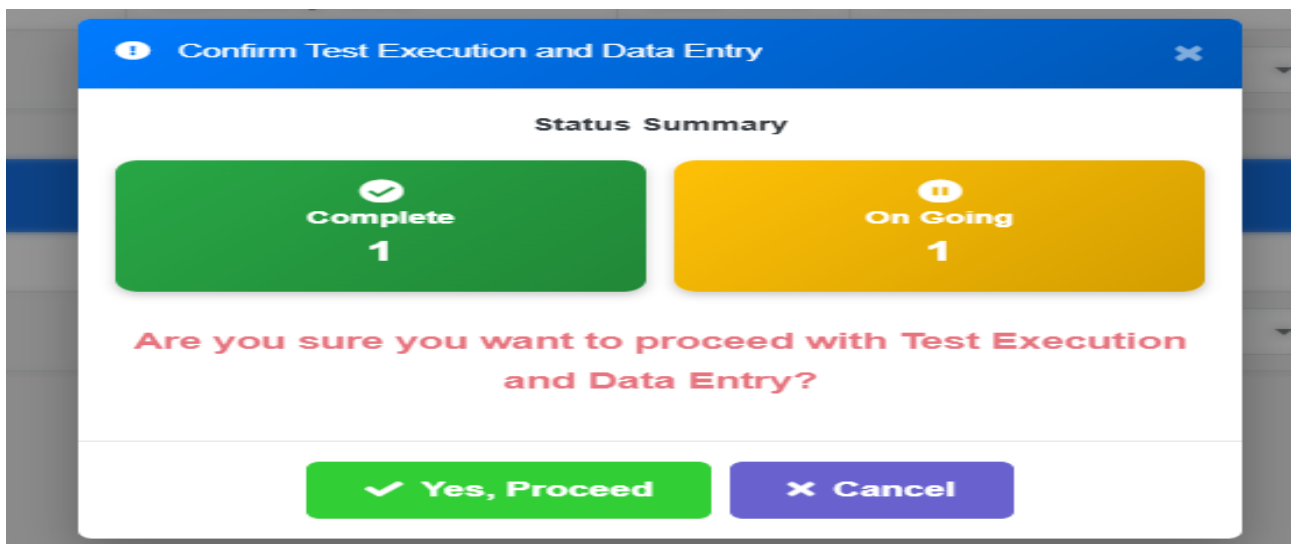

The screenshot shows a 'Confirm Test Execution and Data Entry' modal. At the top, it has a title bar with an information icon and a close button. Below the title bar is a 'Status Summary' section. This section contains two large colored boxes: a green box on the left labeled 'Complete' with a checkmark icon and the number '1', and a yellow box on the right labeled 'On Going' with a pause icon and the number '1'. Below these boxes is a red text prompt asking 'Are you sure you want to proceed with Test Execution and Data Entry?'. At the bottom of the modal are two buttons: a green 'Yes, Proceed' button and a purple 'Cancel' button.

Confirm Test Execution and Data Entry

Status Summary

Complete
1

On Going
1

Are you sure you want to proceed with Test Execution and Data Entry?

✓ Yes, Proceed **✕ Cancel**

Data Review Module

The Data Review module allows users to perform a critical review of the results against test parameters and can also validate the result to move them forward or reject them with remarks.

Data Review > Transactional > Data Review

All Fields

Job Order No: MMCQA-JN-0126-0002 Samples: 2

No.	Date Submitted	Laboratory Code	Sample Name	Laboratory Tests	Lead Time	Status	Lab Result	Remarks	Review Verification	Logs
1	January 02, 2026	12-0126-0003	ARGENTINE CORN	CRUDE FAT	3	COMPLETE	2231	No Remarks	Select	View Logs
2	January 02, 2026	12-0126-0003	ARGENTINE CORN	NIR PROXIMATE	2	COMPLETE	321312	No Remarks	Select	View Logs

Click the dropdown and select the review verifications then click save, confirmation modal will appear click “Yes, proceed” to move forward to Result Verification Module.

Confirm Data Review

Status Summary

Verified
1

Re-Analysis
1

Are you sure you want to proceed with Data Review?

Yes, Proceed

Cancel

Note: *If Re-analysis is selected the results will proceed to Test Execution and Data Entry.*

Result Verification Module

Final Checkpoint before report released. Users with access can perform this module.

Result Verification Home > Transactional > Result Verification

All Fields Search here

Job Order No: MMCQA-JN-0126-0002 Samples: 1

						Test Execution		Data Review	Result Verification		
No.	Date Submitted	Laboratory Code	Sample Name	Laboratory Tests	Lead Time	Status	Lab Result	Remarks	Review Verificat...	Test Result	Remarks
1	January 02, 2026	12-0126-0003	ARGENTINE CORN	NIR PROXIMATE	2	COMPLETE	321312	No Remarks	VERIFIED	Select	

Click the dropdown test results, *if selected is disapproved the result will proceed to Test Execution and Data Entry and Remarks will be required.* Click “Save” then confirmation modal will appear click “Yes, proceed” to move format to Report Releasing Module.

Confirm Result Verification

Status Summary

 Approved 1	 Disapproved 0
--	---

Are you sure you want to proceed with Result Verification?

✓ Yes, Proceed **✗ Cancel**

Report Releasing Module

Allow users to Report Released the approved results and can also download here the Certificate of Analysis.

Report Release [Home](#) > [Transaccional](#) > [Report Release](#)

All Fields

Job Order No: MMCQA-JN-0126-0002
Client: JAMES DIMACULANGAN
Nutritionist: MR. ALFREDO BAUTISTA

Samples: 1 [Select All](#)

No.	Date Submitted	Laboratory Code	Sample Name	Laboratory T...	Result	Test Execution	Data Review	Result Ver...	With COA re...	For Releasing	Reference No.	Generate COA	Logs
1	January 02, 2026	12-0126-0003	ARGENTINE CORN	NIR PROXIMATE	321312	COMPLETE	VERIFIED	APPROVED	Yes	☐	-	Generate COA	View Logs

Click the checkbox of the laboratory code you want to release then click “Save” confirmation modal will appear click “Yes, proceed” to release the laboratory code.

Confirm Releasing

Status Summary

✓
Number of Releasing
1

Are you sure you want to proceed with Releasing?

✓ **Yes, Proceed** ✕ **Cancel**

Here you can see we are able to released the laboratory code and the system automatically generate the reference number of that COA then click “Generate COA” button to download it.

Report Release
Transactional > Report Release

All Fields
Search here

Job Order No: MMCQA-JN-0126-0002
Client: JAMES DIMACULANGAN
Nutritionist: MR. ALFREDO BAUTISTA

Samples: 1
Select All

No.	Date Submitted	Laboratory Code	Sample Name	Laboratory T...	Result	Test Execution	Data Review	Result Ver...	With COA re...	For Releasing	Reference No.	Generate COA	Logs
1	January 02, 2026	12-0126-0003	ARGENTINE CORN	NIR PROXIMATE	321312	COMPLETE	VERIFIED	APPROVED	Yes		26-MMCQA-00001		View Logs

This is Certificate Analysis Sample.



MAMBATANGAN MILLING CORPORATION QUALITY ASSURANCE SATELLITE
LABORATORY
MAMBATANGAN, MANOLO FORTICH, BUKIDNON 8703

CERTIFICATE OF ANALYSIS

Reference No. 26-MMCQA-00001

Customer : BOUNTY PLUS INC.
Address : SAMPLE ADDRESS
Tel./Fax No : 09221311231
Laboratory Code : 12-0126-0003
Sample Name : ARGENTINE CORN
Date Received : Jan 02, 2026
Date Analyzed : Jan 02, 2026
Date Reported : Jan 02, 2026

TEST PARAMETER	RESULTS	TEST METHOD
CRUDE FIBER, %	321312	NIR PREDICTION (BRUKER)

Remarks:
Results of analysis as per sample submitted. Samples will be kept only for a month from the date received.
*SEP – Standard Error of Prediction

REFERENCE/S:
TANGO FOR REFLECTION AND MEASUREMENT USER MANUAL

Analyzed by:

JAMES DIMACULANGAN
Registered Chemical Technician
System Admin
License No.: 102021
Valid Until: November 19, 2025

Certified True and Correct:

RUDEL LIM
Registered Chemist
QA Laboratory Supervising Chemist
License No.: 0013476
Valid Until: September 20, 2028

Signed for the Company by:

MAYETTE DOMENCIL
Registered Chemist
QA Laboratory Manager
License No.: 008870
Valid Until: May 01, 2026

Best Practices

- Verify test values before saving
- Log out after use
- Do not share login credentials

Support & Contact Information

For assistance contact your system administrator.