

# BLIMS

## Bounty Laboratory Information Management System

---

### User Manual

#### Introduction

BLIMS (Bounty Laboratory Information Management System) is a web-based system designed to help laboratory personnel record, manage, and monitor laboratory test results efficiently. This manual is intended for non-technical users responsible for recording laboratory sample tests.

#### System Objective.

To modernize request handling, testing workflows, and data security by replacing manual, paper-based processes with a streamlined, digital platform.

#### User Roles & Access

Laboratory users can log in, view assigned samples, record and update test results, and review test history. Access depends on permissions assigned by the system administrator.

#### System Requirements

- Desktop or laptop computer
- Internet connection (with VPN)
- Google Chrome (recommended)

#### Logging In to BLIMS


1. Open your web browser.
2. Enter the BLIMS system URL.
3. Enter your username and password.
4. Click Login.



The image shows a login form for the Bounty Laboratory Info Management System (BLIMS). The form is centered on a white background with a dark blue border. At the top, the title "Bounty Laboratory Info Management" is in bold black text, with "BLIMS" in blue text below it. Below the title are two input fields: "Email or Employee No." and "Password". The "Password" field has a toggle icon (an eye) to its right. A large blue "Login" button is centered below the input fields. At the bottom, the text "BLIMS © 2025" is displayed.

## Registration Overview

The purpose of the Registration module is to enable users to create, view, and maintain sample registration records within the system. To access go to side bar menu and click “Registration”

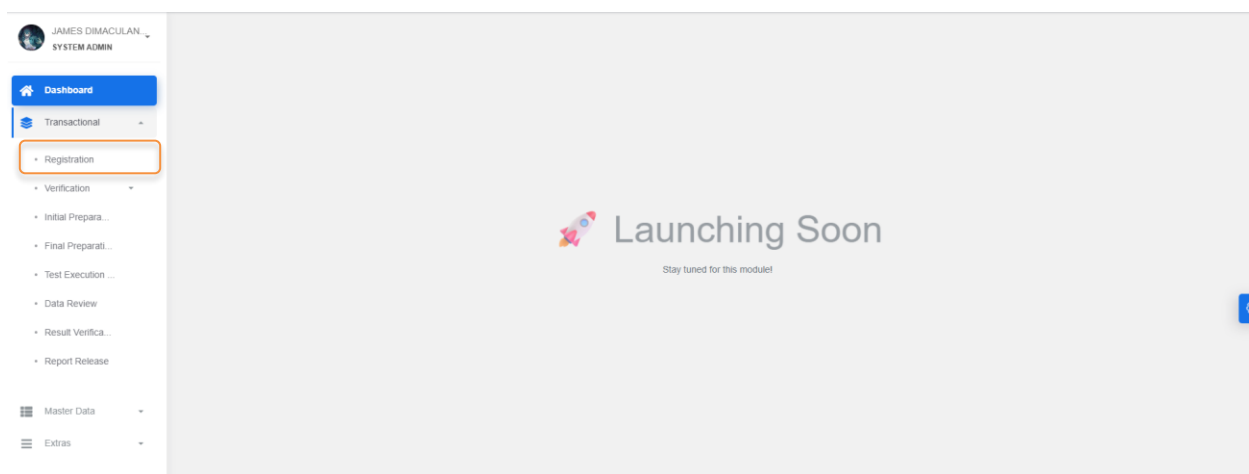

Registration > Transactional > Registration

[Add Registration](#)

All Fields

Job Order No: PFMQA-JN-0126-0001 [Add Detail](#) Samples 2

| No. | Laboratory Code | Sample Name | Production ... | Shipment Supplier   | Plate / ... | Batch / ... | Type o...               | Labora...        | Test P...        | Test Status | Reason              | COA R... | Remarks    | Action |
|-----|-----------------|-------------|----------------|---------------------|-------------|-------------|-------------------------|------------------|------------------|-------------|---------------------|----------|------------|--------|
| 1   | 29-0126-0001    | ADISODIUM   | 2026-01-02     | A & S INTERNATIONAL | 1221        | 2001        | FEEDS AND RAW MATERIALS | MOISTURE CONTENT | MOISTURE, %      | ONGOING     | No Reason Indicated | No       | No Remarks |        |
| 2   | 29-0126-0002    | ADISODIUM   | 2026-01-02     | A & S INTERNATIONAL | ABB 5515    | 2001        | FEEDS AND RAW MATERIALS | CRUDE PROTEIN    | CRUDE PROTEIN, % | ONGOING     | No Reason Indicated | Yes      | No Remarks |        |

Click “Add Registration” to create a new sample registration then you will be redirect to creating registration module.

Registration > Transactional > Registration

Registration

Internal Feedmill \*  
Select

Commercial Feedmill  
Select Feed

Address  
Enter address

External Client

Contact Number \*  
0917 123 4567

Nutritionist \*  
Select

Laboratory Location \*  
Select Location

Delivery Type  
☒ Courier ☐ Internal Handling

Way Bill Number  
Enter Way Bill Number

Submitted By  
JAMES DIMACULANGAN

Date Submitted  
2026-01-02

Remarks  
Enter any remarks

Attach File  
[Choose File](#) No file chosen

Sample Details \*

| Sample Name | Type of Sample | Test Code | Production/Delivery D... | Shipment/Supplier | Plate/Van Number | Batch/Lot Number | Type of Lead Time | COA                      | Action |
|-------------|----------------|-----------|--------------------------|-------------------|------------------|------------------|-------------------|--------------------------|--------|
| Sample      | Type of Sample | Test Code | dd/mm/yyyy               | Supplier          | Plate/Van        | Batch/Lot Number | Select Lead Time  | <input type="checkbox"/> |        |

[Save](#) [Cancel](#)

Here just fill up the all the dropdown details required for all fields.

Note. To enabled the Sample Details, you need to select the feedmills, nutritionist and Laboratory Location.


Sample Details \*

| Sample Name | Type of Sample | Test Code | Production/Delivery D... | Shipment/Supplier | Plate/Van Number | Batch/Lot Number | Type of Lead Time | COA                      | Action |
|-------------|----------------|-----------|--------------------------|-------------------|------------------|------------------|-------------------|--------------------------|--------|
| Sample      | Type of Sample | Test Code | dd/mm/yyyy               | Supplier          | Plate/Van        | Batch/Lot Number | Select Lead Time  | <input type="checkbox"/> |        |

[Save](#) [Cancel](#)

Each field in the sample details is a dropdown fields, just select all the fields and also, we have an action button where we can add other samples for this Registration. Click “Save” to Bounty Laboratory Management System (BLIMS) – User Guide

submit. Once save Success Message will appear.

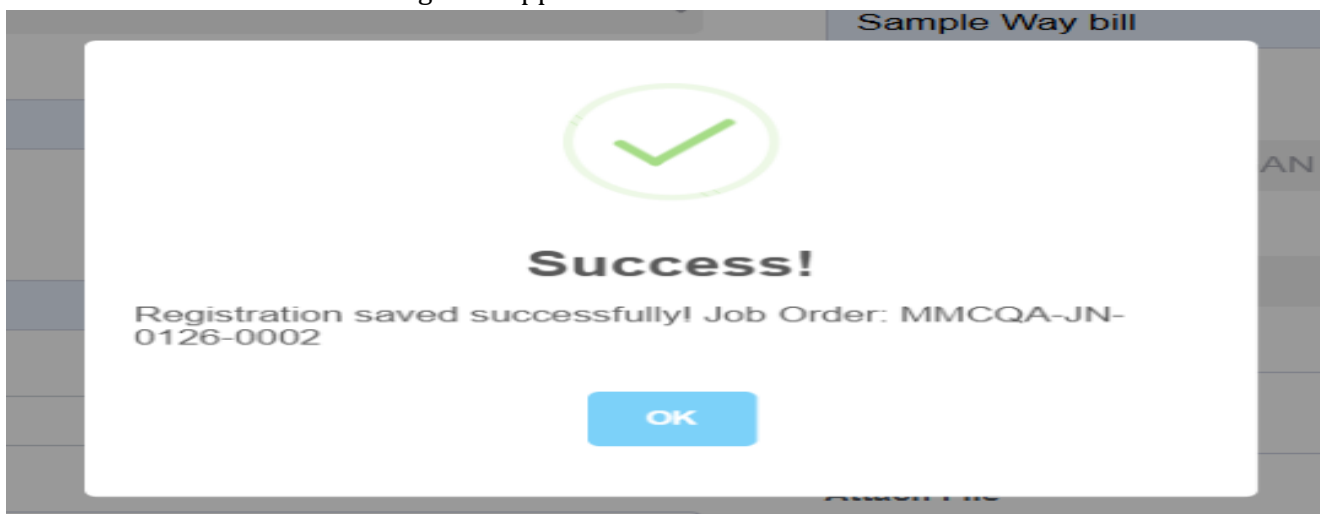

Then go back to View of Created Registration. Here we can edit a registered samples(*only if you have an access to do so*).

| Job Order No: MMCQA-JN-0126-0002 |                 |                |                |                     |             |             |                         |           |              |             |                     |          |            |                                                                                      |
|----------------------------------|-----------------|----------------|----------------|---------------------|-------------|-------------|-------------------------|-----------|--------------|-------------|---------------------|----------|------------|--------------------------------------------------------------------------------------|
|                                  |                 |                |                |                     |             |             |                         |           |              |             |                     |          |            |                                                                                      |
| No.                              | Laboratory Code | Sample Name    | Production ... | Shipment Supplier   | Plate / ... | Batch / ... | Type o...               | Labora... | Test P...    | Test Status | Reason              | COA R... | Remarks    | Action                                                                               |
| 1                                | 12-0126-0003    | ARGENTINE CORN | 2026-01-02     | A & S INTERNATIONAL | AJA 9142    | B#C2523 -07 | FEEDS AND RAW MATERIALS | CRUDE FAT | CRUDE FAT, % | ONGOING     | No Reason Indicated | Yes      | No Remarks |  |

Click “Add Detail” if you want to add a new samples to existing Job Order then modal will appear fill up all the details then Click “Save” if done.

Add Sample Detail

Sample Name

Select

Type of Sample

Select

Test Code

Select

Production/Delivery Date

dd/mm/yyyy

Shipment/Supplier

Select

Plate / Van Number

Select

Batch / Lot Number

Select

Type of Lead Time

COA Required:

☐

Save

Cancel

If you want to append a new sample with sample lab code under the same Job Order, Click the Action button “replicate icon” then modal will appear and it will displayed the same details you want to replicate.

Job Order No: MMCGA-JN-0126-0002

+ Add Detail

Samples: 1

| No. | Laboratory Code | Sample Name    | Production ... | Shipment Supplier   | Plate / ... | Batch / ... | Type o...                | Labora... | Test P...    | Test Status | Reason              | COA R... | Remarks    | Action      |
|-----|-----------------|----------------|----------------|---------------------|-------------|-------------|--------------------------|-----------|--------------|-------------|---------------------|----------|------------|-------------|
| 1   | 12-0126-0003    | ARGENTINE CORN | 2026-01-02     | A & S INTERNATIONAL | AJA 9142    | B#C2523 -07 | FEEDS AND RAW MATERI ALS | CRUDE FAT | CRUDE FAT, % | ONGOING     | No Reason Indicated | Yes      | No Remarks | <div></div> |

Replicate Sample Details: 12-0126-0003

Sample Name

ARGENTINE CORN

Type of Sample

Select

Test Code

Select

Production/Delivery Date

02/01/2026

Shipment/Supplier

A & S INTERNATIONAL

Plate / Van Number

AJA 9142

Batch / Lot Number

B#C2523-07

Type of Lead Time

COA Required:

☒

Save

Cancel

Then fill up all the required fields then Click “Save”, then modal for confirmation will appear, click “Yes, Proceed”.

Replicate Sample Details: 12-0126-0003

Sample Name

ARGENTINE CORN

Type of Sample

LIMESTONE

Test Code

NIR-CF

Production/Delivery Date

02/01/2026

Shipment/Supplier

A & S INTERNA

Plate / Van Number

AJA 9142

Batch / Lot Number

B#C2523-07

Type of Lead Time

Regular

COA Required:

☒

Save

Cancel

Replicating this Sample

Are you sure you want to Replicate this Samples?


✓ Yes, Proceed

✕ Cancel


## Verification Module

Allows users to view and update the verification status of registered samples, categorized as Ongoing, Cancelled, or Failed. Users can assign verification results and provide reasons or


remarks when required.

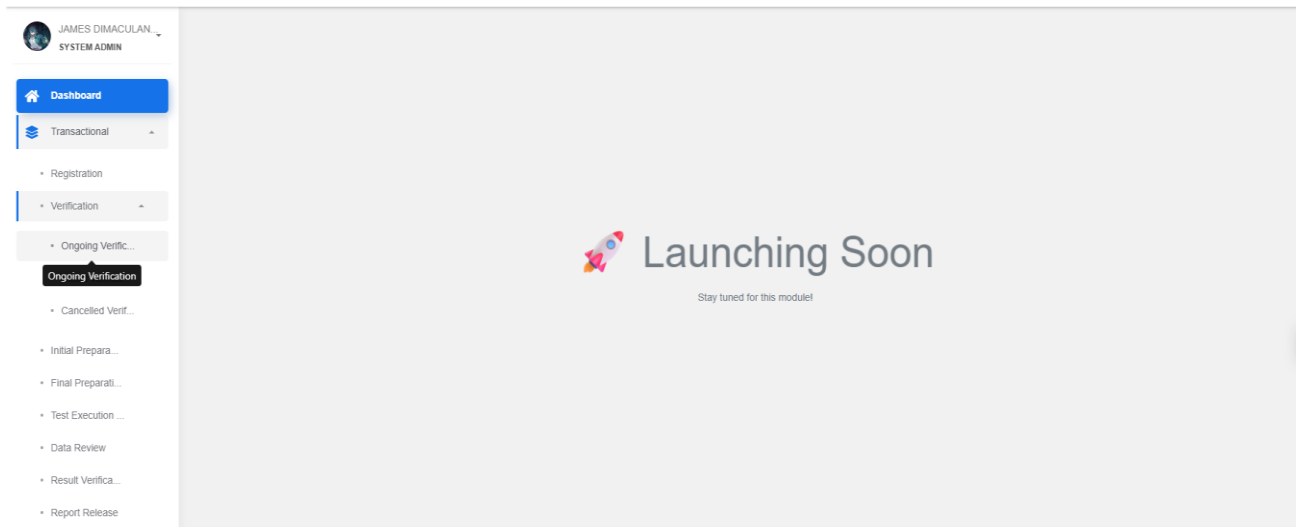

Click Ongoing Verification Module.

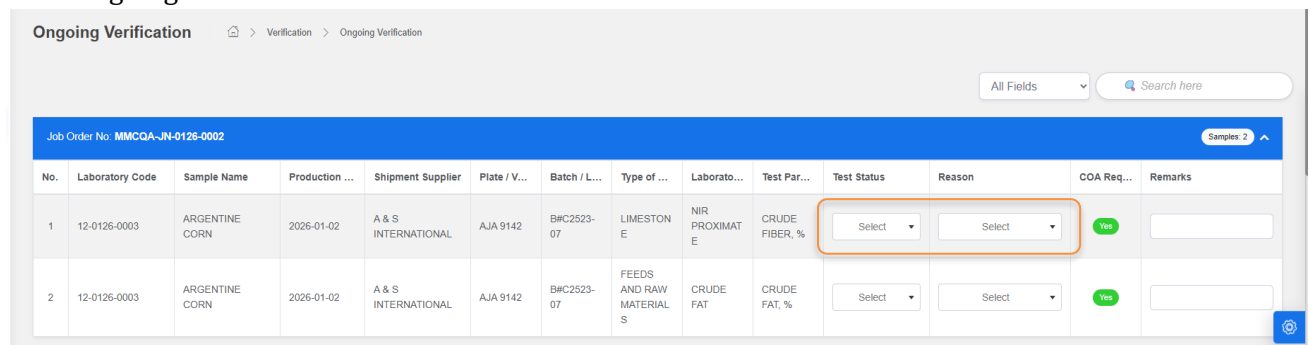

Select test status and reasons for this. *If On hold and failed the selected test status, reasons and remarks field are required. It will also notify the User who created the Registration via email.* Then Click “Save” then Confirmation modal will appear.

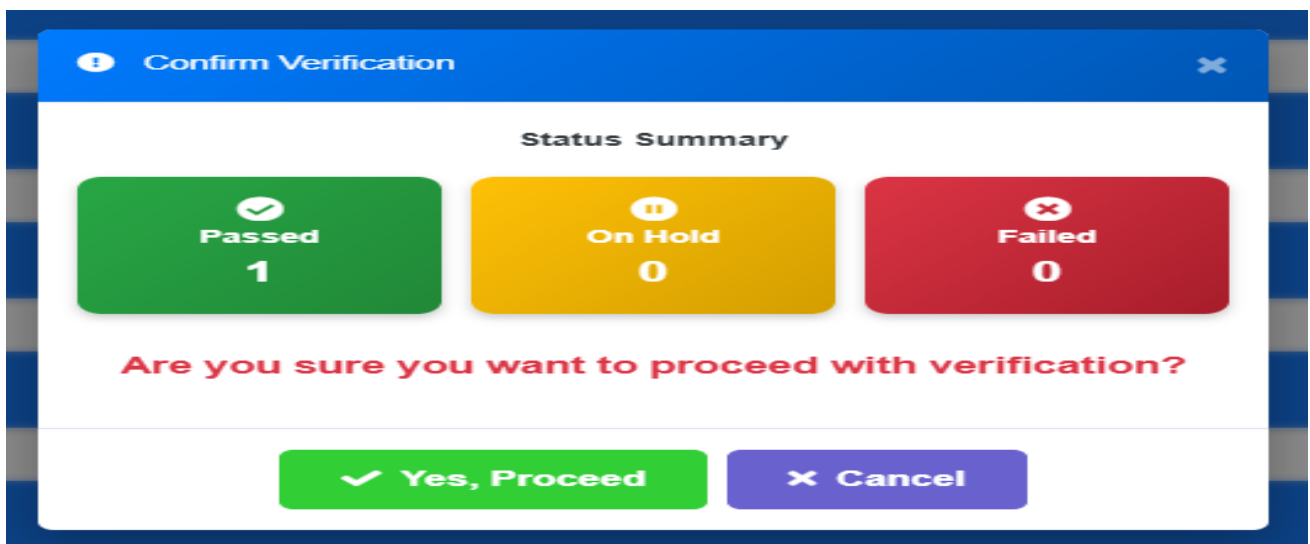



Click “Yes, proceed”.

**Failed Verification** [Home](#) > [Verification](#) > [Failed Verification](#)

All Fields

Job Order No: PFMQA-JN-1125-0054 Samples: 1

| No. | Laboratory Code | Sample Name    | Production ... | Shipment Supplier | Plate / V... | Batch / L... | Type of ... | Laborato...   | Test Par... | Test Status | Reason                        | COA Req... | Remarks |
|-----|-----------------|----------------|----------------|-------------------|--------------|--------------|-------------|---------------|-------------|-------------|-------------------------------|------------|---------|
| 1   | 16-1125-0068    | ARGENTINE SOYA | 2025-11-26     | AAFCO             | ABB 5515     | BWC2523-08   | MOLASSE S   | NIR PROXIMATE | ASH, %      | FAILED      | DAMAGED OR IMPROPER PACKAGING | No         | failed  |

*Note: If samples are not edit or change test status after 3 business days it will automatically Cancelled by our system.*

**Cancelled Verification** [Home](#) > [Verification](#) > [Cancelled Verification](#)

All Fields

Job Order No: MMCQA-JN-1125-0089 Samples: 3

| No. | Laboratory Code | Sample Name | Production ... | Shipment Supplier   | Plate / V... | Batch / L... | Type of ... | Laborato...   | Test Par...    | Test Status | Reason                                                              | COA Req... | Remarks                                        |
|-----|-----------------|-------------|----------------|---------------------|--------------|--------------|-------------|---------------|----------------|-------------|---------------------------------------------------------------------|------------|------------------------------------------------|
| 1   | 16-1125-0085    | ADISODIUM   | 2025-11-27     | A & S INTERNATIONAL | 1221         | 2001         | LIMESTONE   | NIR PROXIMATE | CRUDE FIBER, % | CANCELLED   | DAMAGED OR IMPROPER PACKAGING                                       | No         | Auto-cancelled by system after 3 business days |
| 2   | 16-1125-0089    | ADISODIUM   | 2025-11-27     | A & S INTERNATIONAL | 1221         | BWL2323-02   | LIMESTONE   | NIR PROXIMATE | CRUDE FIBER, % | CANCELLED   | INCORRECT SAMPLE WEIGHT                                             | No         | Auto-cancelled by system after 3 business days |
| 3   | 16-1225-0002    | ADISODIUM   | 2025-11-27     | ACX DIAMOND         | 1221         | BWC2523-07   | OIL         | CA(EDTA)      | CALCIUM, %     | CANCELLED   | INSUFFICIENT RESOURCES (E.G., CHEMICAL REAGENTS, CONSUMABLES, ETC.) | No         | Auto-cancelled by system after 3 business days |

## Initial Preparation Module

The Initial Preparation module allows users to manage and edit the pre-analytical process before moving samples to Final Preparation.


**Initial Preparation** [Home](#) > [Transactional](#) > [Initial Preparation](#)

All Fields

Job Order No: MMCQA-JN-0126-0002 Samples: 1

| No. | Laboratory Code | Sample Name    | Production ... | Shipment Supplier   | Plate / Van Number | Batch / Lot Number | Type of Sample | Laboratory Tests | Test Parameter | Pre Analytical Process |
|-----|-----------------|----------------|----------------|---------------------|--------------------|--------------------|----------------|------------------|----------------|------------------------|
| 1   | 12-0126-0003    | ARGENTINE CORN | 2026-01-02     | A & S INTERNATIONAL | AJA 9142           | BWC2523-07         | LIMESTONE      | NIR PROXIMATE    | CRUDE FIBER, % | Select                 |

Select Pre-analytical Process for this Sample then Click “Save” then confirmation modal will appear. Click “Yes, proceed”

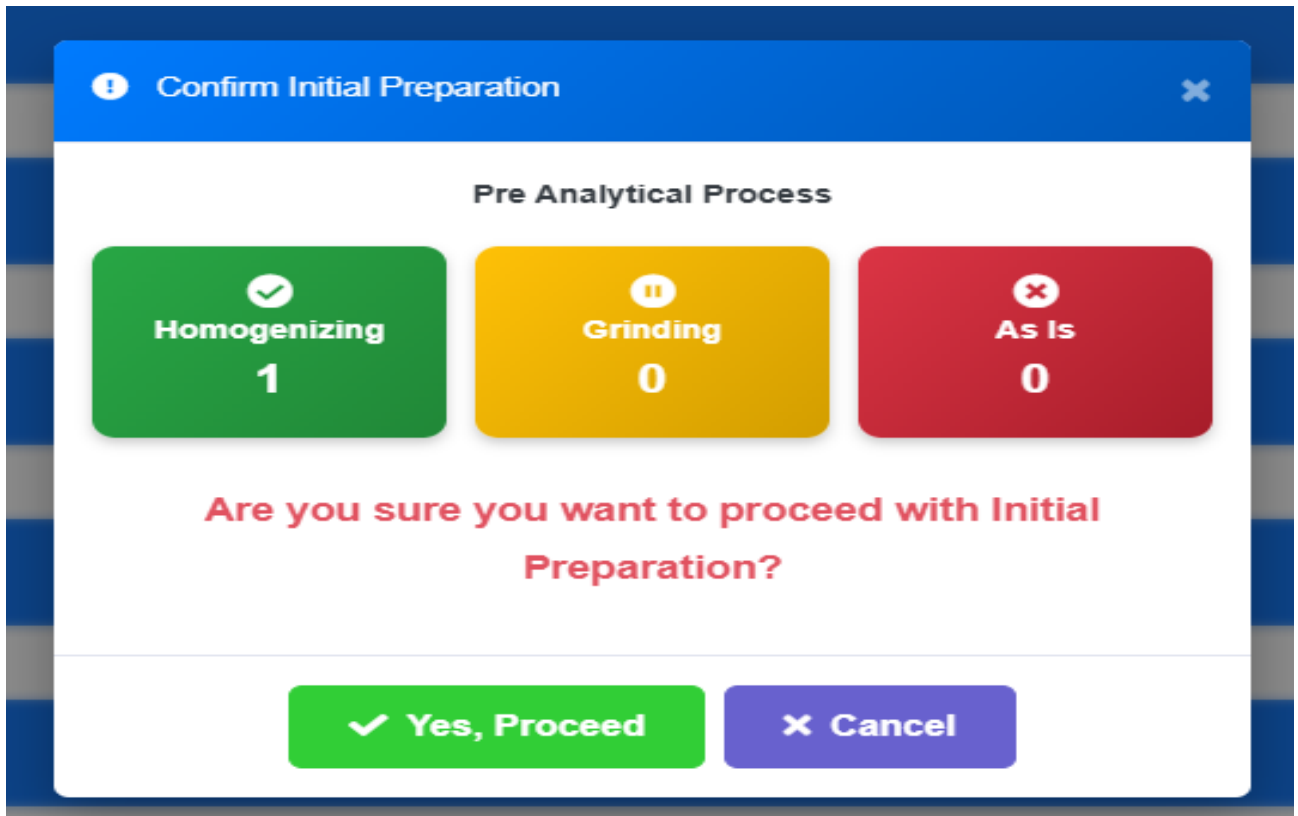

## Final Preparation Module

The Final Preparation module allows users to manage and edit the preparation Verification before moving samples to Test Execution and Data Entry.


Final Preparation Home > Transactional > Initial Preparation > Final Preparation

All Fields Search here

Job Order No: MMCQA-JN-0126-0002 Samples: 1

| No. | Laboratory Code | Sample Name    | Production ... | Shipment Suppl...   | Plate / Van Num... | Batch / Lot Nu... | Type of Sample | Laboratory Tests | Test Parameter | Pre Analytical A... | Preparation Ver... | Remarks |
|-----|-----------------|----------------|----------------|---------------------|--------------------|-------------------|----------------|------------------|----------------|---------------------|--------------------|---------|
| 1   | 12-0126-0003    | ARGENTINE CORN | 2026-01-02     | A & S INTERNATIONAL | AJA 9142           | BAC2523-07        | LIMESTONE      | NIR PROXIMATE    | CRUDE FIBER, % | HOMOGENIZING        | Select             |         |

Select in the dropdown field. *If failed remarks are required and it will notify the user via email and the samples will proceed to Initial Preparation.* Click "Save" then formation modal will appear.

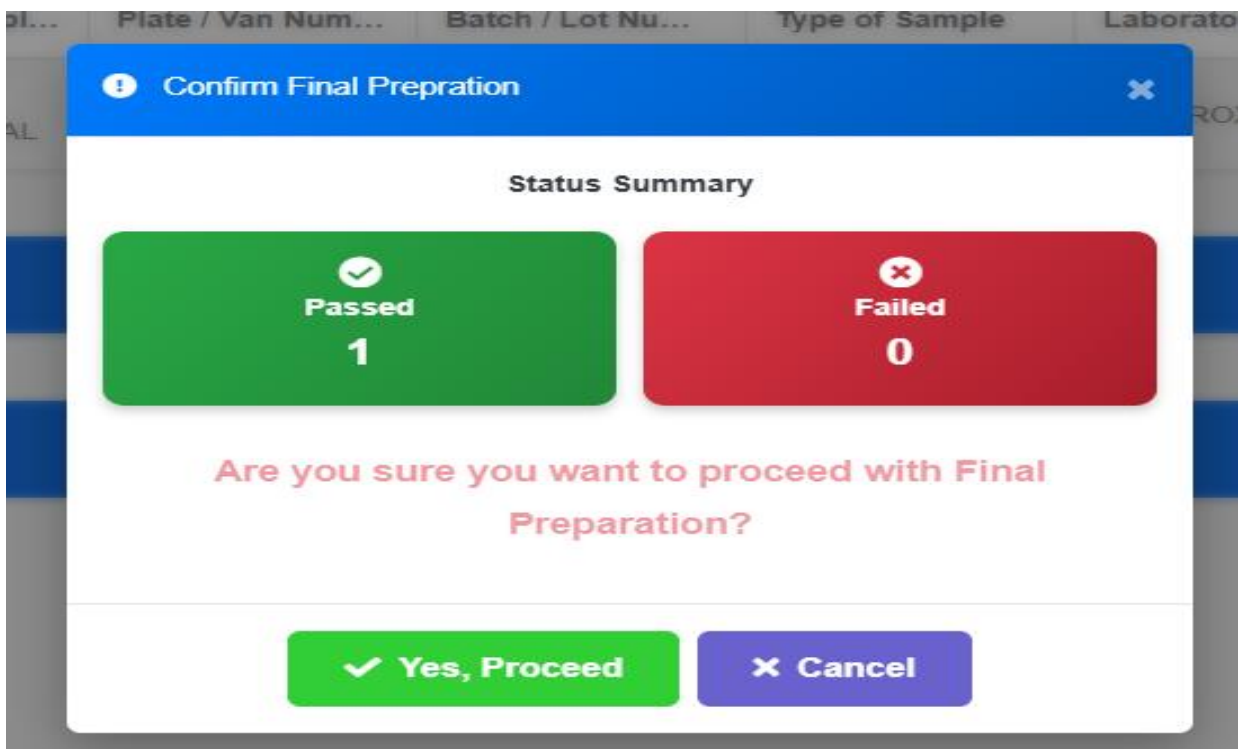

**Confirm Final Preparation**

**Status Summary**


|                                                                                                                |                                                                                                                  |
|----------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------|
| <br><b>Passed</b><br><b>1</b> | <br><b>Failed</b><br><b>0</b> |
|----------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------|

**Are you sure you want to proceed with Final Preparation?**

Click “Yes, proceed” to move the samples to Test Execution and Data Entry.

## Test Execution and Data Entry Module

The Test Execution and Data Entry module allows users to manage and edit the status also provide the Lab Result and Remarks; it also displayed here the Lead time.

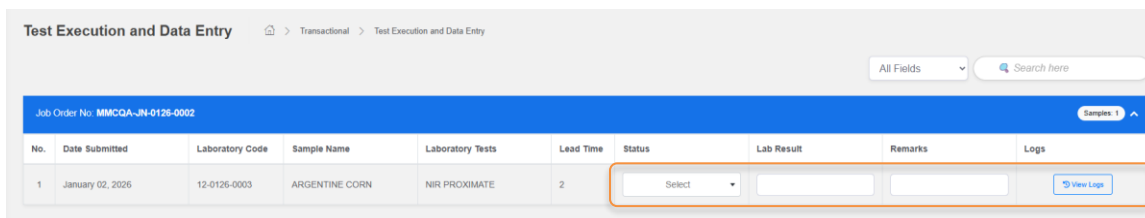

**Test Execution and Data Entry**

Job Order No: MMCQA-JN-0126-0002


| No. | Date Submitted   | Laboratory Code | Sample Name    | Laboratory Tests | Lead Time | Status                              | Lab Result           | Remarks              | Logs                                     |
|-----|------------------|-----------------|----------------|------------------|-----------|-------------------------------------|----------------------|----------------------|------------------------------------------|
| 1   | January 02, 2026 | 12-0126-0003    | ARGENTINE CORN | NIR PROXIMATE    | 2         | <input type="text" value="Select"/> | <input type="text"/> | <input type="text"/> | <input type="button" value="View Logs"/> |

Select for dropdowns.

If the selected is Completed in the status Lab result is required.

If the Sample Beyond Lead time it will notify Manager, Supervising Chemist and Analyst via email and Remarks is required.


Click “View logs” to view samples history or record.

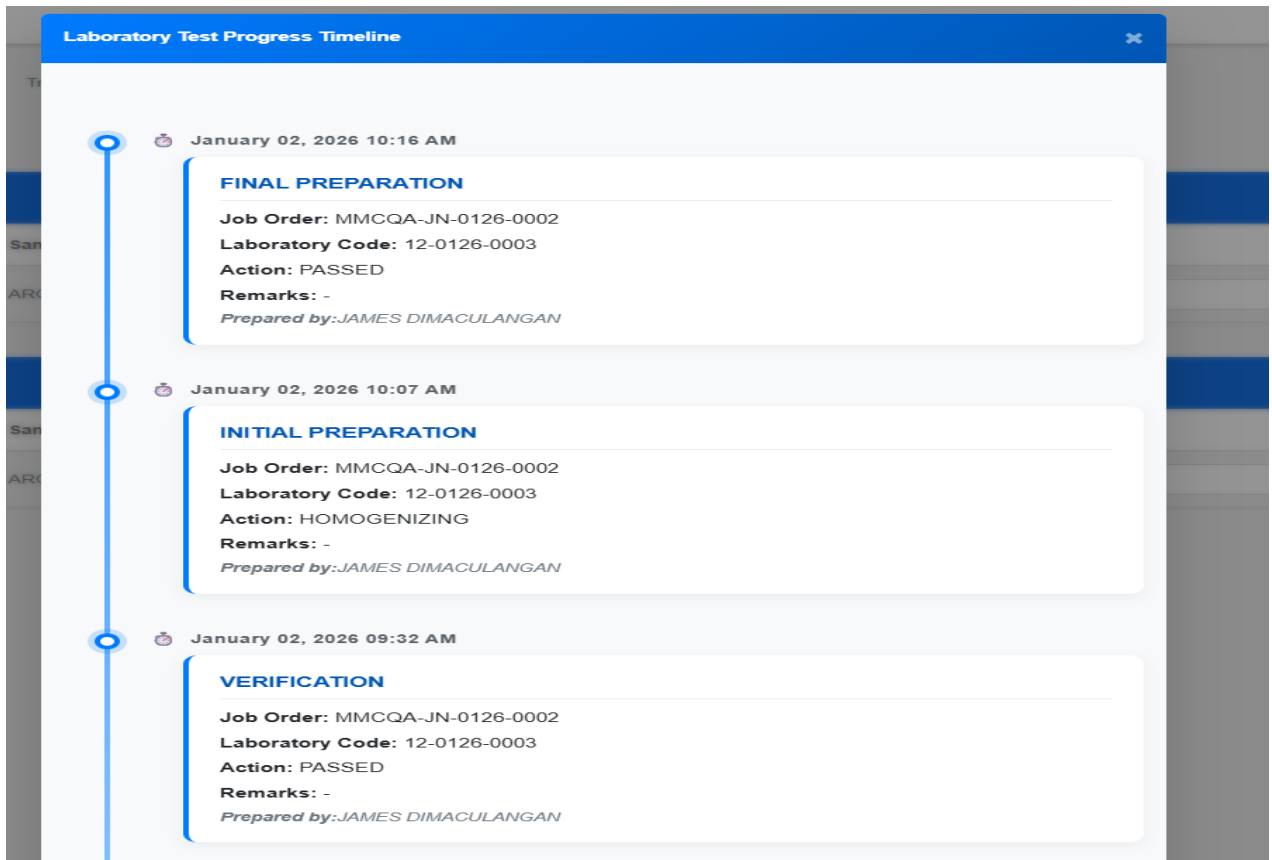

**Laboratory Test Progress Timeline**

- January 02, 2026 10:16 AM**  
**FINAL PREPARATION**  
**Job Order:** MMCQA-JN-0126-0002  
**Laboratory Code:** 12-0126-0003  
**Action:** PASSED  
**Remarks:** -  
*Prepared by: JAMES DIMACULANGAN*
- January 02, 2026 10:07 AM**  
**INITIAL PREPARATION**  
**Job Order:** MMCQA-JN-0126-0002  
**Laboratory Code:** 12-0126-0003  
**Action:** HOMOGENIZING  
**Remarks:** -  
*Prepared by: JAMES DIMACULANGAN*
- January 02, 2026 09:32 AM**  
**VERIFICATION**  
**Job Order:** MMCQA-JN-0126-0002  
**Laboratory Code:** 12-0126-0003  
**Action:** PASSED  
**Remarks:** -  
*Prepared by: JAMES DIMACULANGAN*

Click “Save”, then confirmation modal will appear click “Yes, Proceed” to move the samples for Data Review Module.

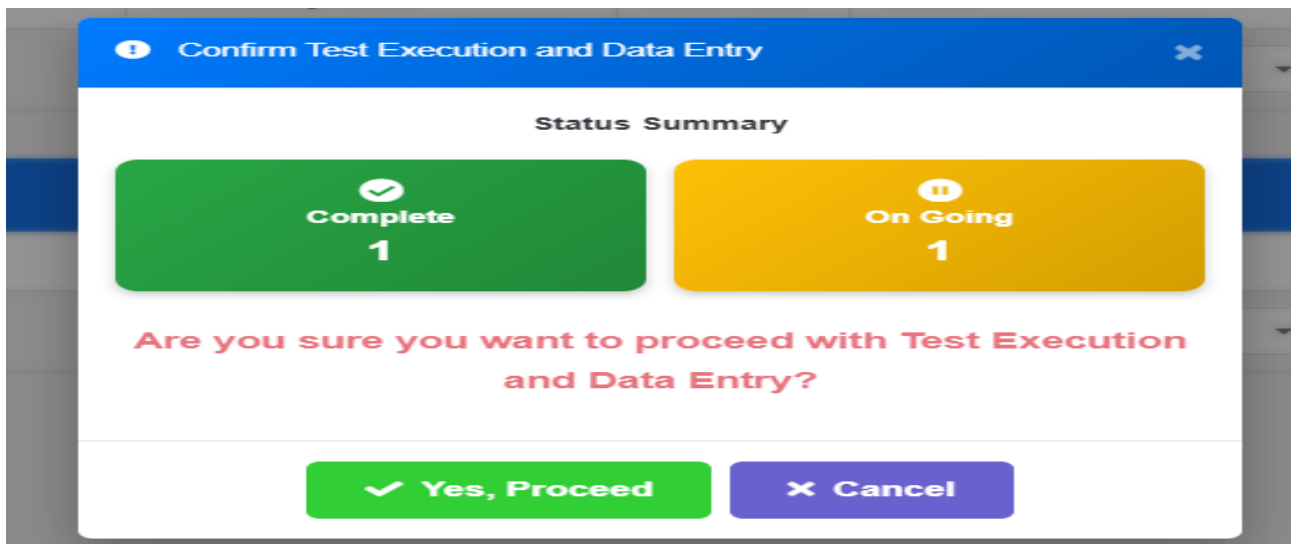

**Confirm Test Execution and Data Entry**

**Status Summary**

|                                                                                                                    |                                                                                                                      |
|--------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------|
| <br><b>Complete</b><br><b>1</b> | <br><b>On Going</b><br><b>1</b> |
|--------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------|

**Are you sure you want to proceed with Test Execution and Data Entry?**

## Data Review Module

The Data Review module allows users to perform a critical review of the results against test parameters and can also validate the result to move them forward or reject them with remarks.

**Data Review** [Transactional](#) > [Data Review](#)

All Fields

Job Order No: MMCQA-JN-0126-0002 Samples: 2

| No. | Date Submitted   | Laboratory Code | Sample Name    | Laboratory Tests | Lead Time | Status   | Lab Result | Remarks    | Review Verification | Logs                      |
|-----|------------------|-----------------|----------------|------------------|-----------|----------|------------|------------|---------------------|---------------------------|
| 1   | January 02, 2026 | 12-0126-0003    | ARGENTINE CORN | CRUDE FAT        | 3         | COMPLETE | 2231       | No Remarks | Select              | <a href="#">View Logs</a> |
| 2   | January 02, 2026 | 12-0126-0003    | ARGENTINE CORN | NIR PROXIMATE    | 2         | COMPLETE | 321312     | No Remarks | Select              | <a href="#">View Logs</a> |

Click the dropdown and select the review verifications then click save, confirmation modal will appear click “Yes, proceed” to move forward to Result Verification Module.

**Confirm Data Review**

**Status Summary**

**Verified**  
1

**Re-Analysis**  
1

**Are you sure you want to proceed with Data Review?**

**Yes, Proceed**

**Cancel**

Note: *If Re-analysis is selected the results will proceed to Test Execution and Data Entry.*

## Result Verification Module

Final Checkpoint before report released. Users with access can perform this module.

Result Verification Transaction > Result Verification

Job Order No: MMCQA-JN-0126-0002 Samples 1

| No. | Date Submitted   | Laboratory Code | Sample Name    | Laboratory Tests | Lead Time | Status   | Lab Result | Remarks    | Review Verificat... | Test Result | Remarks |
|-----|------------------|-----------------|----------------|------------------|-----------|----------|------------|------------|---------------------|-------------|---------|
| 1   | January 02, 2026 | 12-0126-0003    | ARGENTINE CORN | NIR PROXIMATE    | 2         | COMPLETE | 321312     | No Remarks | VERIFIED            | Select      |         |

Click the dropdown test results, *if selected is disapproved the result will proceed to Test Execution and Data Entry and Remarks will be required*. Click “Save” then confirmation modal will appear click “Yes, proceed” to move format to Report Releasing Module.

**Confirm Result Verification**

**Status Summary**

**Approved**  
**1**

**Disapproved**  
**0**

**Are you sure you want to proceed with Result Verification?**

**Yes, Proceed**

**Cancel**

## Report Releasing Module

Allow users to Report Released the approved results and can also download here the Certificate of Analysis.


**Report Release** [Home](#) > [Transactional](#) > [Report Release](#)

Job Order No: MMCQA-JN-0126-0002  
Client: JAMES DIMAGULANGAN  
Nutritionist: MR. ALFREDO BAUTISTA

Samples: 1 [Select All](#)

| No. | Date Submitted   | Laboratory Code | Sample Name    | Laboratory T... | Result | Test Execution | Data Review | Result Ver... | With COA re... | For Releasing            | Reference No. | Generate COA             | Logs                      |
|-----|------------------|-----------------|----------------|-----------------|--------|----------------|-------------|---------------|----------------|--------------------------|---------------|--------------------------|---------------------------|
| 1   | January 02, 2026 | 12-0126-0003    | ARGENTINE CORN | NIR PROXIMATE   | 321312 | COMPLETE       | VERIFIED    | APPROVED      | Yes            | <input type="checkbox"/> | -             | <a href="#">Download</a> | <a href="#">View Logs</a> |

Click the checkbox of the laboratory code you want to release then click “Save” confirmation modal will appear click “Yes, proceed” to release the laboratory code.

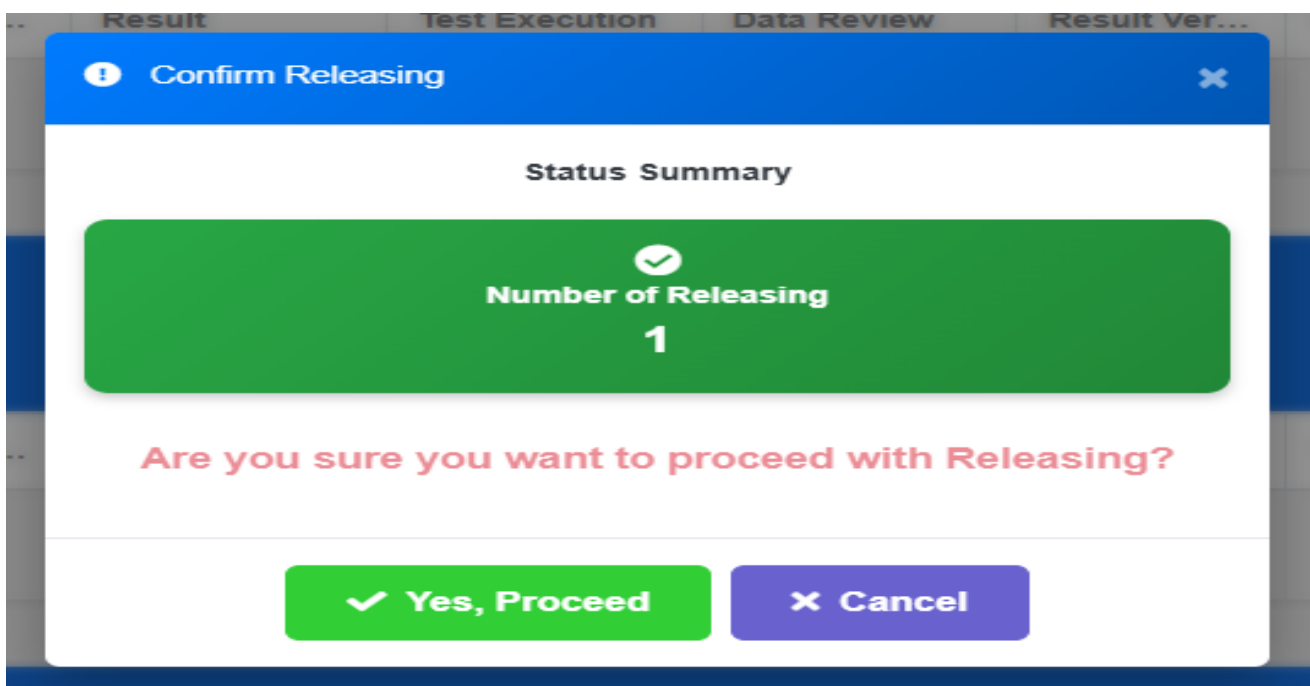

Here you can see we are able to released the laboratory code and the system automatically generate the reference number of that COA then click “Generate COA” button to download it.

Report Release [Home](#) > [Transactional](#) > [Report Release](#)

All Fields

Job Order No: MMCQA-JN-0126-0002  
 Client: JAMES DIMACULANGAN  
 Nutritionist: MR. ALFREDO BAUTISTA

Samples: 1 [Select All](#)

| No. | Date Submitted   | Laboratory Code | Sample Name    | Laboratory T... | Result | Test Execution | Data Review | Result Ver... | With COA re... | For Releasing | Reference No.  | Generate COA | Logs                      |
|-----|------------------|-----------------|----------------|-----------------|--------|----------------|-------------|---------------|----------------|---------------|----------------|--------------|---------------------------|
| 1   | January 02, 2026 | 12-0126-0003    | ARGENTINE CORN | NIR PROXIMATE   | 321312 | COMPLETE       | VERIFIED    | APPROVED      | Yes            |               | 26-MMCQA-00001 |              | <a href="#">View Logs</a> |

This is Certificate Analysis Sample.



**MAMBATANGAN MILLING CORPORATION QUALITY ASSURANCE SATELLITE  
 LABORATORY  
 MAMBATANGAN, MANOLO FORTICH, BUKIDNON 8703**

### CERTIFICATE OF ANALYSIS

Reference No. 26-MMCQA-00001

Customer : **BOUNTY PLUS INC.**  
 Address : SAMPLE ADDRESS  
 Tel./Fax No : 09221311231  
 Laboratory Code : **12-0126-0003**  
 Sample Name : **ARGENTINE CORN**  
 Date Received : Jan 02, 2026  
 Date Analyzed : Jan 02, 2026  
 Date Reported : Jan 02, 2026

| TEST PARAMETER | RESULTS | TEST METHOD             |
|----------------|---------|-------------------------|
| CRUDE FIBER, % | 321312  | NIR PREDICTION (BRUKER) |

**Remarks:**

Results of analysis as per sample submitted. Samples will be kept only for a month from the date received.  
 \*SEP – Standard Error of Prediction

**REFERENCE/S:**

TANGO FOR REFLECTION AND MEASUREMENT USER MANUAL

Analyzed by:



**JAMES DIMACULANGAN**  
 Registered Chemical Technician  
 System Admin  
 License No.: 102021  
 Valid Until: November 19, 2025

Certified True and Correct:



**RUDEL LIM**  
 Registered Chemist  
 QA Laboratory Supervising Chemist  
 License No.: 0013476  
 Valid Until: September 20, 2028

Signed for the Company by:

**MAYETTE DOMENCIL**  
 Registered Chemist  
 QA Laboratory Manager  
 License No.: 008870  
 Valid Until: May 01, 2026



## Best Practices

- Verify test values before saving
- Log out after use
- Do not share login credentials

## Frequently Asked Questions (FAQs)

### 1. The system is running slowly. What should I do?

If you are experiencing system slowdown:

- Contact your **System Administrator** or **IT Technical personnel** assigned to your area.
- Ensure your internet connection is stable.

### 2. The system cannot be accessed. What could be the reason?

If you cannot access the system, possible reasons include:

- **Slow or unstable internet connection**, as BLIMS is an online system.
- **Incorrect URL** entered in the browser. Make sure the correct system URL is used (e.g., *site domain URL*).
- **Server downtime**, where the system server is temporarily unavailable.

If the issue persists, contact your **System Administrator** or **IT Support**.

### 3. Some modules or transactions are missing. Why?

If you notice that certain modules or transactions are not visible:

- You may **not have access permission** to those modules.
- Access rights depend on your **assigned user role**.

To resolve this:

- Verify with your **System Administrator** if you should have access.
- Request the necessary **permission or role update** if required.

### 4. The screen layout overlaps or does not display properly when zoomed in. What should I do?

Bounty Laboratory Management System (BLIMS) – User Guide

BLIMS uses tables with multiple columns. When the browser zoom level is set too high, the screen layout may overlap or display incorrectly.

**Recommendations:**

- Set your browser zoom level to **100% or 90%**.
- Use **horizontal scrolling** instead of increasing the zoom level.
- **Maximize the browser window** to improve screen visibility.
- Access the system using a **desktop or laptop computer** for best results.

**Note:** Increasing the zoom level may affect how tables and columns are displayed in the system.

## **5. Who should I contact for system support?**

For any system-related concerns, such as access issues, errors, or performance problems:

- Contact your **System Administrator**
- Or reach out to the **IT Technical personnel** assigned to your area

## **Support & Contact Information**

For assistance contact your system administrator.